

VAGINAL HIFU

HIGH INTENSE FOCUSED ULTRASOUND



F
A
S
T
S
E
R
V
I
C
E
A
N
D
B
E
S
T
V
A
L
U
E



VAGINAL TIGHTENING HIFU
10000 shots



How the HIFU vaginal tightening works?


Introduction

The Vagina Tightening HIFU System uses a noninvasive ultrasonic focusing technique to directly focus on the mucosal lamia and muscle layer. Using ultrasonic waves as the energy source and taking advantage of its penetration and focality, the system will send out ultrasonic energy focusing in the the lamina propria and muscle fiber layer in a predetermined depth. A higher intensity of ultrasonic region, called focus region, is formed. In 0.1 second, the temperature of the region can reach to above 65 °C , so the collagen is reorganized and the normal tissue outside the focal region is undamaged. Therefore, the desired depth layer can obtain the ideal effect of collagen contraction, reorganization and regeneration. Ultimately, the mysterious effect of vagina tightening is achieved.



Technique Creates Beauty

The newest technology for Vaginal Rejuvenation

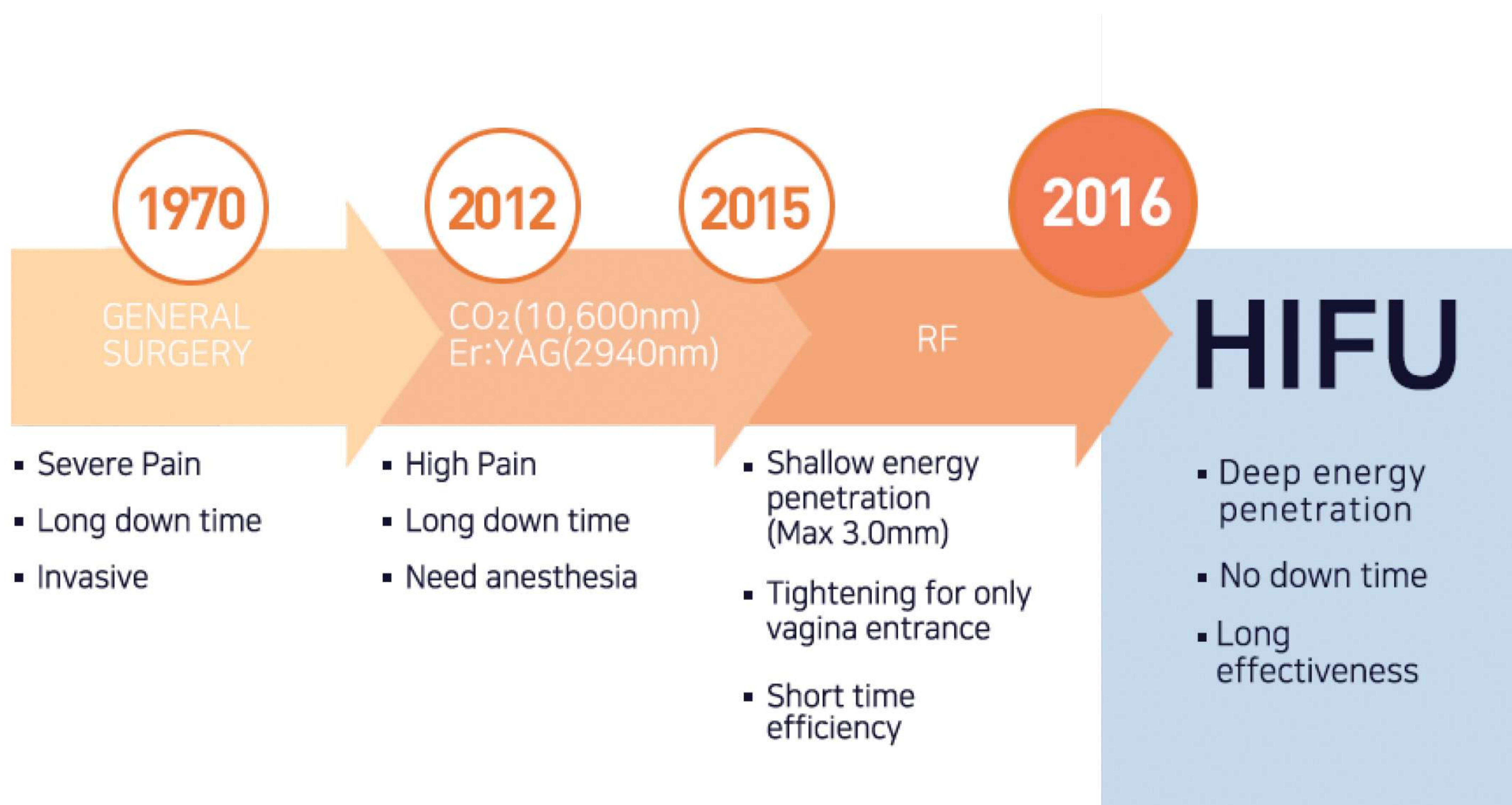


Benefits of Vaginal HIFU

- Increased firming and tightening effect
- Improved sexual Life
- Increased Vaginal lubrication
- Improved arousal sensitivity by cell-regeneration
- Reduces Vaginal Infections
- No dryness
- 360° automatic control of rotation angle

Mechanism of New Collagen Regeneration

In 3 weeks, collagen remodeling process is completed. Also, it will last for 1 year through the result of collagen neogenesis and remodeling

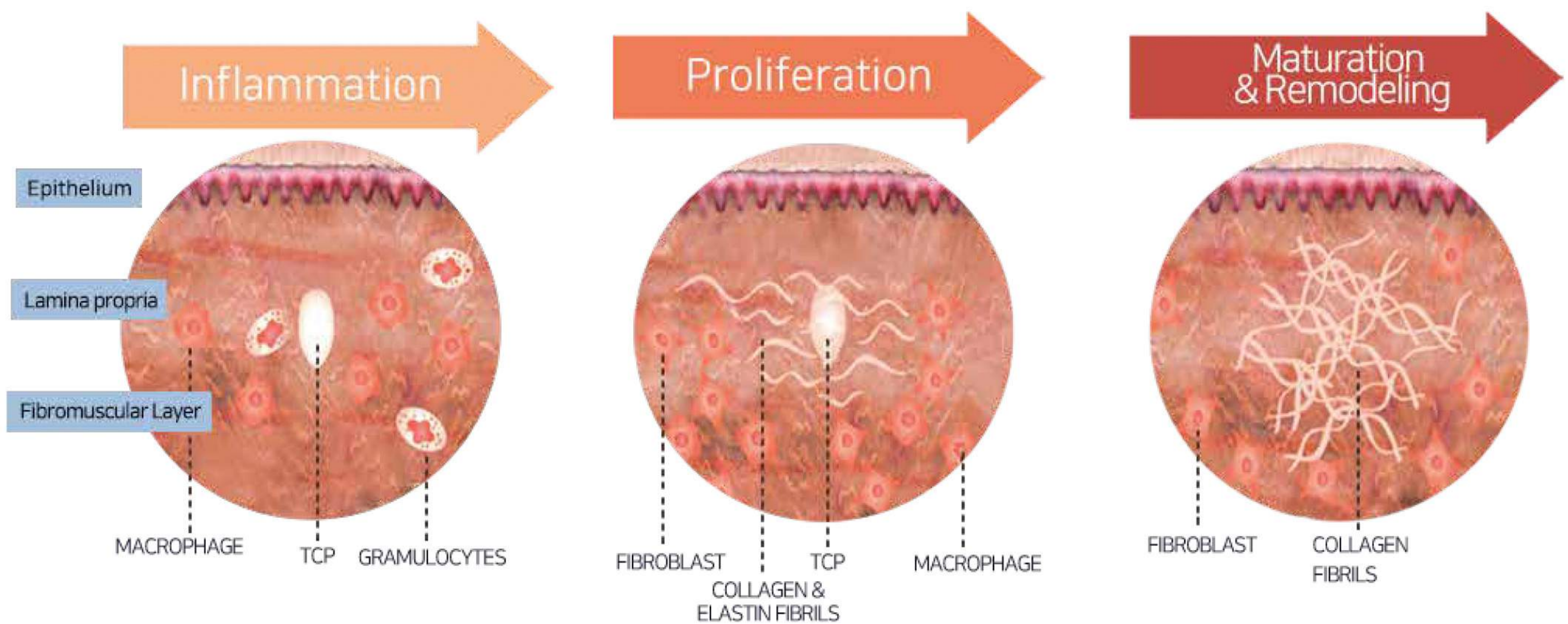

Indications for Vaginal HIFU

- Post menopause
- Vaginal Tightening
- Eliminating dryness
- Recurrent Infections
- SUI: Stress Urinary Incontinence
- Vaginal Axis
- Rugae
- Vaginal lubrication
- Mucosal integrity
- Collagen Regeneration

Technique Creates Beauty

2 Cartridges For Different Depth Treatment



CARTRIDGE
M3.0
7.0MHz, 3.0mm


CARTRIDGE
D4.5
4.0MHz, 4.5mm



- ✓ **4.5mm heads:** Produce the energy directly to SMAS, make it heat coagulation, make the SMAS tighten and lifting. Let the muscles improve the muscle structure from deep to shallow, make better to help the muscle layer restore elasticity and tighten;
- ✓ **3.0mm heads:** Produce energy to subnormal tissue accelerate cell activity and make collagen regrowth and restructuring, increase mucous elasticity and make vaginal muscle tighten.

Technique Creates Beauty

DETAIL DISPLAY

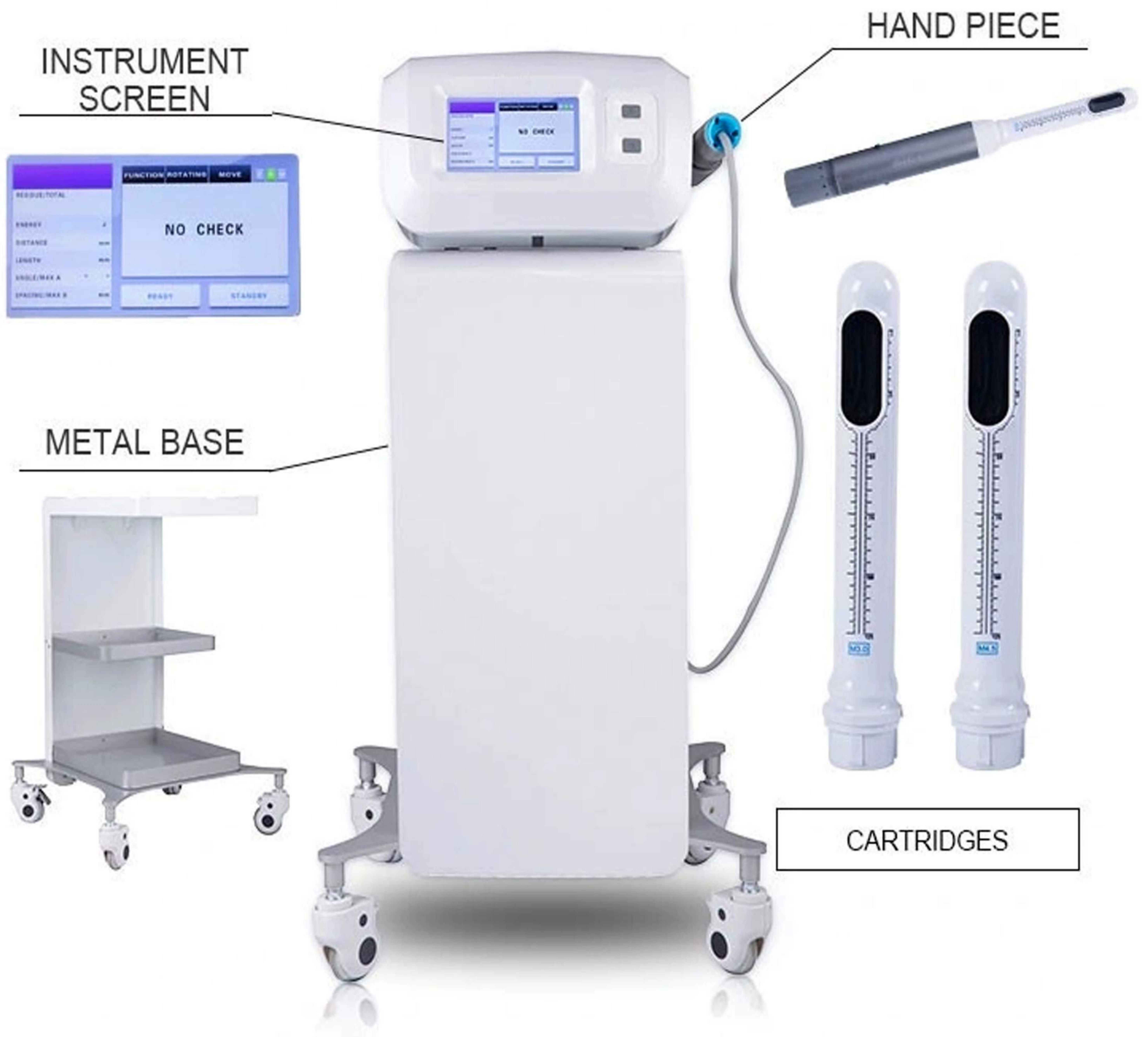

Technique Creates Beauty

■ Technique Specification

Deep	Accurate reach 3.0mm or 4.5mm
Temperature	65°C~70°C
Effect	Collagen Regeneration
Maintain time	1 year
Damage	Non-invasive
Pain	Low
Epidermis protection	Energy into the dermis layer of skin without injury
Treatment	Interval of 6 months